

Spot-tag® antibody [28A5]

Product code: 28a5


Properties

Description	Mouse monoclonal antibody to Spot-tag®.
Target / Specificity	Spot-tag® (PDRVRAVSHWSS)
Species Specificity	N/a
GenBank Accession Number	N/a
GeneID (NCBI)	N/a
Product Type	Primary antibody
Format	Monoclonal mouse antibody
Source	Mouse
Isotype / Subclass	Mouse IgG1
Clonality	Monoclonal
Clone	28a5
RRID	AB_2892245
Immunogen	Spot-tagged fusion protein
Conjugate	Unconjugated
Purification Method	Protein A purification
Form	Liquid
Concentration	1000 µg/mL
Tested Applications	WB, ELISA
Positive Controls	WB: Spot-tagged GFP (N- or C-terminal) added to HEK-293T cells
Cited Applications	
Recommended Dilutions	WB: 1:1,000 - 1:5,000 ELISA: 1:5,000 - 1:10,000
Storage Buffer	PBS, 50% glycerol
Storage Conditions	Aliquot upon receipt and store at -20°C/-4°F. Avoid freeze-thaw cycles.
Stability	Stable for 1 year at -20°C/-4°F after shipment.
Shipment	Shipped on ice
Size	20 µL; 100 µL


Selected validation data

Western Blot

- Blocking: 5% skimmed milk in PBS-T (PBS + 0.075% Tween-20) for 10-30 min at +4°C/+40°F or room temperature
- Washing: PBS-T
- Incubation: dilution of Spot-tag® antibody [28A5] 1:1,000 - 1:5,000 in 5% milk in PBS-T, 1h at room temperature or overnight at +4°C/+40°F
- *Note:* The optimal dilution depends on the application and should be determined by the user. A titration from range from 1:1,000 up to 1:5,000 is recommended.
- *Note:* Image acquisition time may have to be optimized.

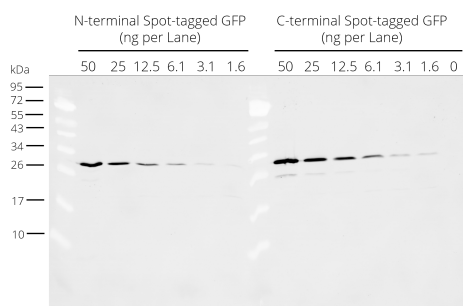

Western blot analysis of N- and C-terminal Spot-tagged GFP added to HEK-293T cell lysate. Detection with Spot-tag® antibody [28A5] (28a5, ChromoTek) 1:5,000 and Nano-Secondary® alpaca anti-mouse IgG1, recombinant VHH, Alexa Fluor® 488 [CTK0103, CTK0104] (sms1AF488-1, ChromoTek) 1:5,000. Western blot membrane was incubated simultaneously with the primary antibody and Nano-Secondary® (one-step staining).

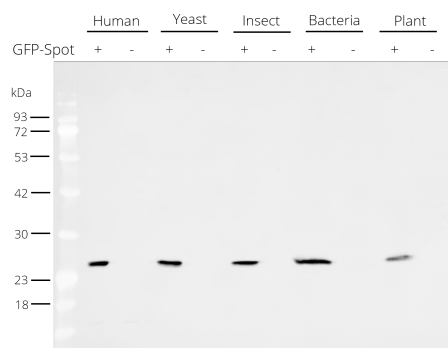

Western blot analysis of different cell lysates +/- 50 ng C-terminal Spot-tagged GFP: human (HEK-293T), yeast (*S.cerevisiae*), insect (High Five™), bacteria (*E.coli*), and plant (*A.thaliana*). Detection with Spot-tag® antibody [28A5] (28a5, ChromoTek) 1:5,000 and Nano-Secondary® alpaca anti-mouse IgG1, recombinant VHH, Alexa Fluor® 488 [CTK0103, CTK0104] (sms1AF488-1, ChromoTek) 1:5,000. Western blot membrane was incubated simultaneously with the primary antibody and Nano-Secondary® (one-step staining).

Spot-tag® antibody [28A5]

Product code: 28a5



Product overview and related products

Spot-tag® toolbox	Product code
Spot-Trap® Agarose	eta-10; -20; -100
Spot-Trap® Agarose Kit	etak-20
Spot-Trap® Magnetic Agarose	etma-10; -20; -100
Spot-Trap® Magnetic Agarose Kit	etmak-20
Spot-Trap® Magnetic Particles M-270	etd-10; -20; -100
Spot-Trap® Magnetic Particles M-270 Kit	etdk-20
iST Spot-Trap® Kit for IP/MS	etak-iST-8
Binding Control Agarose	bab-20
Binding Control Magnetic Agarose	bmab-20
Spin columns	sct-10; sct-20; sct-50
Spot-tag® antibody [28a5]	28a5-20; -100
Spot-Label® Alexa Fluor® 488	ebAF488-10; -50
Spot-Label® Alexa Fluor® 568	ebAF568-10; -50
Spot-Label® Alexa Fluor® 647	ebAF647-10; -50
Spot-Label® ATTO488	eba488-10; -50
Spot-Label® ATTO594	eba594-10; -50
Spot VHH, recombinant binding protein (bivalent)	etb-250
Spot-Cap®	eca-2
Spot-peptide	ep-1; -10
Spot-Cap® and peptide	eca-ep
Spot Vectors for cloning: pSpot1 vector, E. coli, Spot-tag N-term., Kan., high expression pSpot2 vector, E. coli, Spot-tag C-term., Kan., high expression pSpot3 vector, E. coli, Spot-tag C-term., Amp., low expression pSpot4 vector, E. coli, Spot-tag N-term., Amp., low expression pSpot5 vector, S. cerevisiae, Spot-tag N-term., Leu, CEN, low expression pSpot6 vector, S. cerevisiae, Spot-tag C-term., Leu, CEN, low expression pSpot7 vector, S. cerevisiae, Spot-tag N-term., Leu, 2μ, high expression pSpot8 vector, S. cerevisiae, Spot-tag C-term., Leu, 2μ, high expression	ev-1 ev-2 ev-3 ev-4 ev-5 ev-6 ev-7 ev-8
Spot Vectors - positive controls: pSpot-Tag-Actin vector (plasmid) for expression in mammalian cells pSpot2_GFPSpot-Tag vector (plasmid) for expression in E. coli pSpot8_GFP-Spot-Tag vector (plasmid) for expression in S. cerevisiae	ev-31 ev-32 ev-33

Spot-tag® antibody [28A5]

Product code: 28a5

Nano-Secondary® alpaca anti-mouse IgG1, recombinant VHH, Alexa Fluor® 488 [CTK0103, CTK0104]	sms1AF488-1-10; -100
Nano-Secondary® alpaca anti-mouse IgG1, recombinant VHH, Alexa Fluor® 568 [CTK0103, CTK0104]	sms1AF568-1-10; -100
Nano-Secondary® alpaca anti-mouse IgG1, recombinant VHH, Alexa Fluor® 647 [CTK0103, CTK0104]	sms1AF657-1-10; -100

For product details, information, and ordering visit www.chromotek.com.

Contact

support@chromotek.com

ChromoTek GmbH
Am Klopferspitz 19
82152 Planegg-Martinsried
Germany
phone: +49 89 124 148 80
fax: +49 89 124 148 811

ChromoTek Inc.
62-64 Enter Lane
Islandia, NY 11749
USA
phone: 631 501 1058
fax: 631 501 1060

Disclaimer

Only for research applications, not for diagnostic or therapeutic use!

ChromoTek and GFP-Trap, RFP-Trap, Myc-Trap, Spot-Trap, Spot-Tag, Spot-Label, Spot-Cap, Nano-Secondary, F2H Kit, and Chromobody are registered trademarks of ChromoTek GmbH, part of Proteintech Group. Nano-CaptureLigand and V5-Trap are trademarks of ChromoTek GmbH, part of Proteintech Group. Nanobody is a registered trademark of Ablynx, a Sanofi company. Alexa Fluor is a registered trademark of Life Technologies Corporation, a part of Thermo Fisher Scientific Inc. SNAP-tag is a registered trademark and CLIP-tag is a trademark of New England Biolabs, Inc. Octet is a registered trademark of FortéBio, a Sartorius brand. High Five is a trademark of Invitrogen Corporation, a part of Thermo Fisher Scientific Inc. Other suppliers' products may be trademarks or registered trademarks of the corresponding supplier each. Statements on other suppliers' products are given according to our best knowledge.