

Human CYB5R3 Sandwich ELISA Kit Datasheet

Please read it entirely before use

Catalogue Number: KE00069

Size: 5*96T

Sensitivity: 3.5 pg/mL

Range: 31.25-2000 pg/mL

Usage: For the quantitative detection of human CYB5R3 concentrations in serum, plasma and cell lysate.


This product is for research use only and not for use in human or animal therapeutic or diagnostic.

Table of content	page
1. Background	3
2. Principle	3
3. Required Materials	3
4. Kit Components and Storage	4
5. Safety Notes	4
6. Sample Collection and Storage	4
7. Regent Preparation	5
8. Assay Procedure Summary	6
9. Validation Data	7
9.1 Standard curve	7
9.2 Precision	7
9.3 Recovery	8
9.4 Sample values	8
9.5 Sensitivity	8
9.6 Linearity	9
10. References	9

1. Background

CYB5R3, also named as DIA1 and B5R, belongs to the flavoprotein pyridine nucleotide cytochrome reductase family and has two naturally occurring isoforms. Isoform 1 is anchored to the cytoplasmic side of the endoplasmic reticulum membrane and mitochondrion outer membrane, while isoform 2 is the soluble form found in red blood cells. There are two types of methemoglobinemia CYB5R3-related. In type 1, the defect affects the soluble form of the enzyme, is restricted to red blood cells, and causes well-tolerated methemoglobinemia. In type 2, the defect affects both the soluble and microsomal forms of the enzyme and is thus generalized, affecting red cells, leukocytes and all body tissues. Type 2 methemoglobinemia is associated with mental deficiency and other neurologic symptoms. The ELISA kit is suitable for testing serum, plasma, and cell lysates.

2. Principle

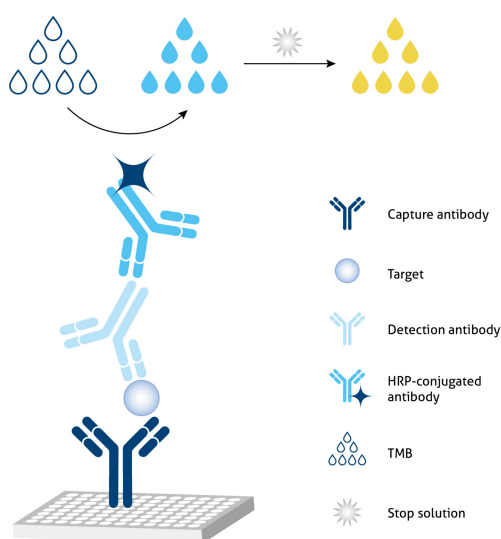

Sandwich ELISA structure (HRP conjugated secondary antibody)

A capture antibody is pre-coated onto the bottom of wells which binds to analyte of interest. A detection antibody also binds to the analyte. Horseradish peroxidase (HRP)-conjugated secondary antibody binds to the detection antibody. TMB acts as the HRP substrate and the solution color will change from colorless to blue. A stop solution containing sulfuric acid turns solution yellow. The color intensity is proportional to the quantity of bound protein which is measurable at 450 nm with the correction wavelength set at 630 nm.

3. Required Materials

- 3.1 A microplate reader capable of measuring absorbance at 450 nm with the correction wavelength set at 630 nm.
- 3.2 Calibrated, adjustable precision pipettes and disposable plastic tips. A manifold multi-channel pipette is recommended for large assays.
- 3.3 Plate washer: automated or manual.
- 3.4 Absorbent paper towels.
- 3.5 Glass or plastic tubes to prepare standard and sample dilutions.
- 3.6 Beakers and graduated cylinders.
- 3.7 Log-log or semi-log graph paper or computer and software for ELISA data analysis. A four-parameter logistic (4-PL) curve-fit is recommended.

4. Kit Components and Storage

Microplate - antibody coated 96-well microplate (8 well × 12 strips)	5 plates	Unopened Kit: Store at 2-8°C for 6 months or -20°C for 12 months. Opened Kit: All reagents stored at 2-8°C for 7 days. Please use a new standard for each assay.
Protein standard - 4000 pg/bottle; lyophilized	10 bottles	
Detection antibody (100×) - 600 μ L/vial*	1 vial	
HRP-conjugated antibody (100×) - 600 μ L/vial*	1 vial	
Sample Diluent PT 1-ec - 150 mL/bottle	1 bottle	
Detection Diluent - 150 mL/bottle	1 bottle	
Wash Buffer Concentrate (20×) - 150 mL/bottle	1 bottle	
Extraction Reagent - 150 mL/bottle	1 bottle	
Tetramethylbenzidine Substrate (TMB) - 60 mL/bottle	1 bottle	
Stop Solution - 60 mL/bottle	1 bottle	
Plate Cover Seals	15 pieces	

* Centrifugation immediately before use

5. Safety Notes

- 5.1 Avoid any skin and eye contact with Stop Solution and TMB. In case of contact, wash thoroughly with water.
- 5.2 Do not use the kit after the expiration date.
- 5.3 Do not mix or substitute reagents or materials from other kit lots or other sources.
- 5.4 Be sure to wear protective equipment such as gloves, masks and goggles during the experiment.
- 5.5 When using an automated plate washer, adding a 30 second soak period following the addition of Wash Buffer to improve assay precision

6. Sample Collection and Storage

- 6.1 Serum: Allow blood samples to clot for 30 minutes, followed by centrifugation for 15 minutes at 1000xg. Clear serum can be assayed immediately or aliquoted and stored at -20°C. Avoid repeated freeze-thaw cycles.
- 6.2 Plasma: Use EDTA, heparin, or citrate as an anticoagulant for plasma collection. Centrifuge for 15 minutes at 1000xg within 30 minutes of collection. The plasma can be assayed immediately or aliquoted and stored at -20°C. Avoid repeated freeze-thaw cycles.
- 6.3 Cell Lysate:
 - 1) Collect cells and wash by centrifuging at 500 x g for 5 minutes before resuspension in pre-cooled PBS buffer. Perform this step three times.
 - 2) Count cells and then discard the supernatant.
 - 3) Add protease inhibitor cocktail to the Extraction Reagent to a final concentration immediately prior to performing cell lysis.
 - 4) Add 1 mL of Extraction reagent (containing protease inhibitor cocktail) Per 1 x 10⁷ cells, Incubate cell suspension on ice for 30 minutes, use ultrasound to treat the samples.
 - 5) Centrifuge cell lysate at 10,000 x g for 5 minutes at 4°C.
 - 6) Measure the concentration of total protein in cell lysate using BCA assay. Where possible, keep samples on ice to avoid protein degradation.

7. Regent Preparation

7.1 Wash Buffer (1X): If crystals have formed in the concentrate, warm to room temperature and mix gently until the crystals have completely dissolved. Add 30 mL of Wash Buffer Concentrate(20X) to 570 mL deionized or distilled water to prepare 1X Wash Buffer.

7.2 Detection Antibody (1X): Dilute 100X Detection Antibody 1:100 using Detection Diluent prior to assay. Suggested 1:100 dilution: 10 μ L 100X Detection Antibody + 990 μ L Detection Diluent (Centrifuge the 100 X Detection Antibody solution for a few seconds prior to use).


7.3 HRP-conjugated antibody (1X): Dilute 100X HRP-conjugated antibody 1:100 using Detection Diluent prior to assay. Suggested 1:100 dilution: 10 μ L 100X HRP-conjugated antibody + 990 μ L Detection Diluent (Centrifuge the 100X HRP-conjugated antibody solution for a few seconds prior to use).

7.4 Sample Dilution: Different samples should be diluted with corresponding Sample Diluent, samples may require further dilution if the readout values are higher than the highest standard OD reading. Variations in sample collection, processing and storage may affect the results of the measurement.

Recommended Dilution for different sample types: 1:2 or 1:4 is recommended for human serum and plasma; 1:4 or 1:8 is recommended for cell lysate.

7.5 Standard Serial Dilution:

Add 2 mL Sample Diluent PT 1-ec in protein standard.

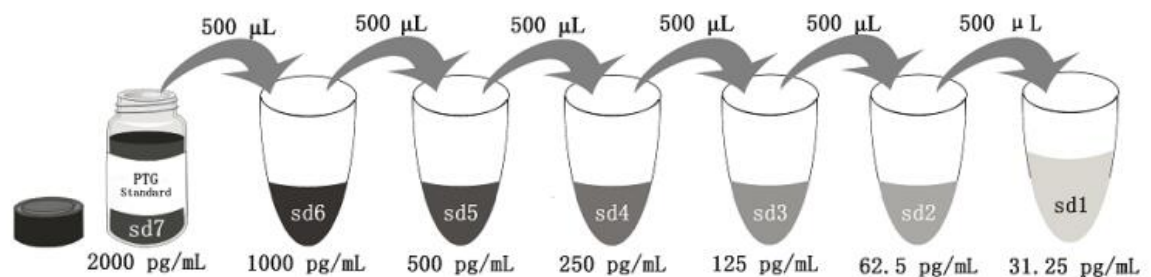

Add # μ L of Standard diluted in the previous step	—	500 μ L	500 μ L	500 μ L	500 μ L	500 μ L	500 μ L
# μ L of Sample Diluent PT 1-ec	2000 μL	500 μ L	500 μ L	500 μ L	500 μ L	500 μ L	500 μ L
	"sd7"	"sd6"	"sd5"	"sd4"	"sd3"	"sd2"	"sd1"

8. Assay Procedure Summary

Bring all reagents to room temperature before use (Detection antibody and HRP-conjugated antibody can be used immediately). To avoid cross-contamination, change pipette tips between additions of each standard level, between sample additions, and between reagent additions. Also, use separate reservoirs for each reagent.

8.1 Take out the required number of microplate strips and return excess strips to the foil pouch containing the drying reagent pack and reseal; store at 4°C immediately. Microplate strips should be used in one week.

8.2 Preset the layout of the microplate, including control group, standard group and sample group, add 100 µL of each standard and sample to the appropriate wells. (Make sure sample addition is uninterrupted and completed within 5 to 10 minutes, it is recommended to assay all standards, controls, and samples in duplicate).

8.3 Seal plate with cover seal, pressing it firmly onto top of microwells. Incubate the plate for 1 hour at 37°C.

8.4 Wash

1) Gently remove the cover seal. Discard the liquid from wells by aspirating or decanting. Remove any residual solution by tapping the plate a few times on fresh paper towels.

2) Wash 4 times with 1X Wash Buffer, using at least 350-400 µL per well. Following the last wash, firmly tap plates on fresh towels 10 times to remove residual Wash Buffer. Avoid getting any towel fibers in the wells or wells drying out completely.

8.5 Add 100 µL of 1X Detection Antibody solution (refer to Reagent Preparation 7.2) to each well. Seal plate with cover seal and incubate for 1 hour at 37°C.

8.6 Repeat wash step in 8.4.

8.7 Add 100 µL of 1X HRP-conjugated antibody solution (refer to Reagent Preparation 7.3) to each well. Seal plate with cover seal and incubate the plate for 40 minutes at 37°C.

8.8 Repeat wash step in 8.4.

8.9 Signal development: Add 100 µL of TMB substrate solution to each well, protected from light. Incubate for 15 to 20 minutes. Substrate Solution should remain colorless until added to the plate.

8.10 Quenching color development: Add 100 µL of Stop Solution to each well in the same order as addition of the TMB substrate. Mix by tapping the side of the plate gently. NB: Avoid skin and eye contact with the Stop solution.

8.11 Read results: Immediately after adding Stop solution read the absorbance on a microplate reader at a wavelength of 450 nm. If possible, perform a double wavelength readout (450 nm and 630 nm).


8.12 Data analysis: Calculate the average of the duplicate readings (OD value) for each standard and sample, and subtract the average of the zero standard absorbance. Construct a standard curve by plotting the mean absorbance for each standard on the y-axis against the concentration on the x-axis, use four-parameter logistic curve-fit (4-PL) analysis to do this. If the samples have been diluted, the OD readout from the standard curve must be multiplied by the dilution factor used.

Step	Reagent	Volume	Incubation	Wash	Notes
1	Standard and Samples	100 µL	60 min	4 times	Cover Wells incubate at 37°C
2	Diluent Antibody Solution	100 µL	60 min	4 times	Cover Wells incubate at 37°C
3	Diluent HRP Solution	100 µL	40 min	4 times	Cover Wells incubate at 37°C
4	TMB Substrate	100 µL	15-20 min	Do not wash	Incubate in the dark at 37°C
5	Stop Solution	100 µL	0 min	Do not wash	-
6	Read plate at 450 nm and 630 nm immediately after adding Stop solution. DO NOT exceed 5 minutes.				

9. Validation Data

9.1 Standard curve

These standard curves are provided for demonstration only. A standard curve should be generated for each set of samples assayed.

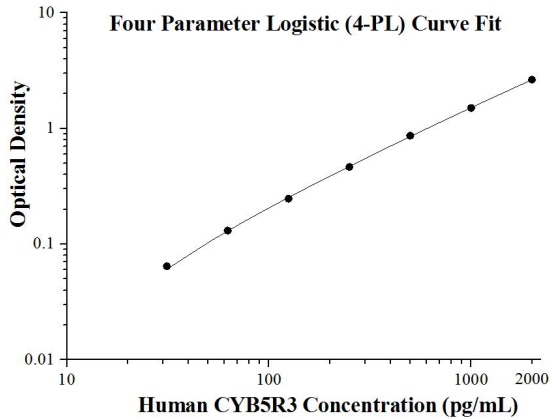

(pg/mL)	O.D	Average	Corrected
0	0.053 0.048	0.0505	-
31.25	0.117 0.113	0.115	0.0645
62.5	0.176 0.187	0.1815	0.131
125	0.294 0.301	0.2975	0.247
250	0.509 0.521	0.515	0.4645
500	0.899 0.935	0.917	0.8665
1000	1.555 1.558	1.5565	1.506
2000	2.403 2.435	2.419	2.3685

9.2 Precision

Intra-assay Precision (Precision within an assay) Three samples of known concentration were tested 20 times on one plate to assess intra-assay precision.

Inter-assay Precision (Precision between assays) Three samples of known concentration were tested in 24 separate assays to assess inter-assay precision.

Intra-assay Precision				
Sample	n	Mean (pg/mL)	SD	CV%
1	20	1,439.3	73.7	5.1
2	20	304.1	18.3	6.0
3	20	75.5	3.3	4.3

Inter-assay Precision				
Sample	n	Mean (pg/mL)	SD	CV%
1	24	1,463.9	76.7	5.2
2	24	313.0	17.4	5.6
3	24	85.1	4.1	4.9

9.3 Recovery

The recovery of human CYB5R3 spiked to three different levels throughout the range of the assay in various matrices was evaluated.

Sample Type		Average% of Expected	Range (%)
Human plasma	1:4	106	100-114
	1:8	107	98-119
Cell lysate	1:50	103	94-105
	1:100	115	112-118

9.4 Sample values

Fourteen serum and plasma samples from healthy volunteers were evaluated for human CYB5R3 in this assay. All samples measured less than the lowest standard, 31.25 pg/mL. No medical histories were available for the donors used in this study.

Cell Lysate:

Sample Type	Concentration (pg/mL)
A549 cell lysate (1 x 10 ⁷ cell)	2,637
MCF-7 cell lysate (1 x 10 ⁷ cell)	6,804

9.5 Sensitivity

The minimum detectable dose of human CYB5R3 is 3.5 pg/mL. This was determined by adding two standard deviations to the concentration corresponding to the mean O.D. of 20 zero standard replicates.

9.6 Linearity

To assess the linearity of the assay, cell lysate samples were diluted with the appropriate **Sample Diluent** to produce samples with values within the dynamic range of the assay. Human plasma samples were spiked with high concentrations of human CYB5R3 in various matrices and diluted with the appropriate **Sample Diluent** to produce samples with values within the dynamic range of the assay. (The samples were initially diluted 1:2)

		Human plasma	Cell Lysate
1:2	Average% of Expected	94	102
	Range (%)	92-96	96-105
1:4	Average% of Expected	101	99
	Range (%)	98-104	97-100
1:8	Average% of Expected	102	114
	Range (%)	98-105	100-125
1:16	Average% of Expected	101	110
	Range (%)	98-105	94-122

10. References

1. Wang Y, Wu YS, et al. A novel mutation in the NADH-cytochrome b5 reductase gene of a Chinese patient with recessive congenital methemoglobinemia. *Blood*. 2000 May 15;95(10):3250-5.
2. Percy MJ, Gillespie MJ, et al. Familial idiopathic methemoglobinemia revisited: original cases reveal 2 novel mutations in NADH-cytochrome b5 reductase. *Blood*. 2002 Nov 15; 100(10): 3447-9
3. Lorenzo FR 5th, Phillips JD, et al. Molecular basis of two novel mutations found in type I methemoglobinemia. *Blood Cells Mol Dis*. 2011 Apr 15;46(4):277-81.