

À des fins de recherche uniquement

Anticorps Polyclonal de lapin anti-Spartin, SPG20

Numéro de catalogue:**CL647-13791** Phare



Informations de base

Numéro de catalogue:	BC047083	Méthode de purification:
CL647-13791		Purification par affinité contre l'antigène
Taille:	Identification du gène (NCBI):	Dilutions recommandées:
100ul , Concentration: 1000 µg/ml by Nanodrop;	23111	IF 1:50-1:500
Hôte:	Nom complet:	Excitation/Emission maxima wavelengths:
Lapin	spastic paraplegia 20 (Troyer syndrome)	654 nm / 674 nm
Isotype:	MW calculé	
IgG	75 kDa, 85 kDa	
Immunogen Catalog Number:	MW observés:	
AG4815	75 kDa, 84 kDa	

Applications

Applications testées:	Contrôles positifs:
FC (Intra), IF	IF : cellules HepG2,
Spécificité de l'espèce:	
Humain, rat, souris	

Informations générales

SPG20 (spastic paraplegia 20) gene encodes a multifunctional Spartn protein. SPG20 protein is highly expressed in adipose tissue and may be implicated in endosomal trafficking and microtubule dynamics. SPG20 is mutated in Troyer syndrome, a hereditary spastic paraplegia defined by degeneration of upper motor neurons, recent study showed that regulation of SPG20 on mitochondrial calcium homeostasis may contribute to the pathophysiology of Troyer syndrome. Hypermethylation of SPG20 promoter, found in colorectal cancer patients, may play a role in cytokinesis arrest in colorectal tumorigenesis.

Stockage


Stockage:
Stocker à -20 °C. Éviter toute exposition à la lumière.
Tampon de stockage:
PBS avec glycérol à 50 %, Proclin300 à 0,05 % et BSA à 0,5 %, pH 7,3.
L'aliquotage n'est pas nécessaire pour le stockage à -20°C

*** Les 20ul contiennent 0,1% de BSA.


For technical support and original validation data for this product please contact:
T: 1(888) 4PTGLAB (1-888-478-4522) (toll free
in USA), or 1(312) 455-8498 (outside USA) E: proteintech@ptglab.com
W: ptglab.com

This product is exclusively available under Proteintech Group brand and is not available to purchase from any other manufacturer.

Données de validation sélectionnées

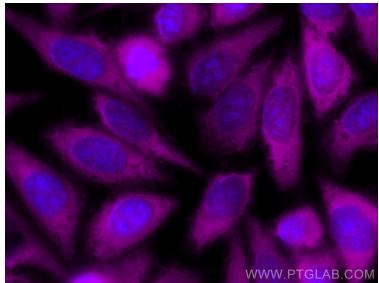

Immunofluorescent analysis of (-20°C Methanol) fixed HepG2 cells using CoraLite® Plus 647 Spartin, SPG20 antibody (CL647-13791) at dilution of 1:200.

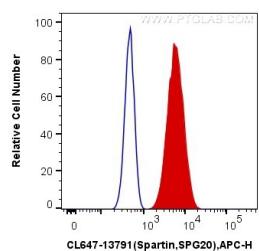

1X10⁶ HepG2 cells were intracellularly stained with 0.4 ug CoraLite® Plus 647 Anti-Human Spartin, SPG20 (CL647-13791) (red), or 0.4 ug Control Antibody. Cells were fixed with 4% PFA and permeabilized with Flow Cytometry Perm Buffer (PF00011-C).