For Research Use Only

CoraLite® Plus 488-conjugated AlaRS Polyclonal antibody

Catalog Number:CL488-17394 Featured Product



Basic Information	Catalog Number: CL488-17394	GenBank Accession Number: BC011451	Purification Method: Antigen affinity purification	
	Size: 100ul, Concentration: 1000 ug/ml by Nanodrop; Source: Rabbit Isotype: IgG Immunogen Catalog Number: AG11294	GeneID (NCBI):	Recommended Dilutions: IF/ICC 1:50-1:500 Excitation/Emission maxima wavelengths: 493 nm / 522 nm	
Applications	Tested Applications: IF/ICC, FC (Intra) Species Specificity: human		Positive Controls: IF/ICC : HepG2 cells,	
Background Information	AARS(alanyl-tRNA synthetase) is also named as AlaRS(alanine tRNA ligase 1, cytoplasmic), renal carcinoma antigen NY-REN-42 and belongs to the class-II aminoacyl-tRNA synthetase family. It can interpret the RNA code and attach specific aminoacids to the tRNAs that contain the cognate trinucleotide anticodons. AARS consists of three domains: the N-terminal catalytic domain, the editing domain and the C-terminal C-Ala domain. The editing domain removes incorrectly charged amino acids, while the C-Ala domain, along with tRNA(Ala), serves as a bridge to cooperatively bring together the editing and aminoacylation centers thus stimulating deacylation of misacylated tRNAs(PMID:19661429).			
Storage	Storage: Store at -20°C. Avoid exposure to ligh Storage Buffer: PBS with 50% Glycerol, 0.05% Proclin Aliquoting is unnecessary for -20°C st	1300, 0.5% BSA, pH 7.3.	ent.	


For technical support and original validation data for this product please contact:T: 1 (888) 4PTGLAB (1-888-478-4522) (toll free
in USA), or 1(312) 455-8498 (outside USA)E: proteintech@ptglab.comW: ptglab.comW: ptglab.com

This product is exclusively available under Proteintech Group brand and is not available to purchase from any other manufacturer.

Selected Validation Data



Immunofluorescent analysis of (-20°C Ethanol) fixed HepG2 cells using CoraLite® Plus 488 AlaRS antibody (CL488-17394) at dilution of 1:200, CL594-Phalloidin (red).

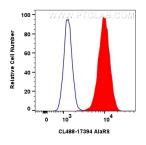

1x10^6 HepG2 cells were intracellularly stained with 0.8 ug Coralite® Plus 488 Alars Polyclonal Antibody (CL488-17394)(red), or 0.8 ug Coralite® Plus 488-conjugated Rabbit IgG control Rabbit PolyAb (CL488-30000) (blue). Cells were fixed with 4% PFA and permeabilized with Flow Cytometry Perm Buffer (PF00011-C).