For Research Use Only

Phantom Dye UV 450 Viability Dye Catalog Number: PD00003



www.ptglab.com

Description Phantom Dye UV 450 is a fixable viability dye which reacts with protein amine

groups in both dead and live cells. Due to the compromised cell membrane in dead cells the dye can bind to many more proteins inside the cell, resulting in much higher fluorescence and easy differentation of Live/Dead cells. Phantom Dye UV 450 is a water soluble UV fluorescent dye excited by the UV laser (355 nm) making it an

excellent addition to your multicolor panel.

Product Information Excitation Laser: UV (355 nm)

Emission (nm): 450

Use By: 6 months from date of receipt

Solution in anhydrous DMSO: 1 l/test in DMSO

Package 1007


Storage -20°C protected from light and moisture.

Application Notes Phantom Dye UV 450 Viability Dye can be used at 1uL/mL of cell suspension.

Concentration required for optimal performance should be determined empirically

by investigator.

Validation Data

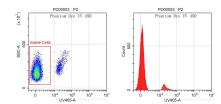

Mouse thymocytes cells were heat killed at 65°C for 10 minutes and then mixed with live mouse thymocytes cells. Cells were then stained with Phantom Dye UV 450 Viability Dye. Viable gate is indicated.